
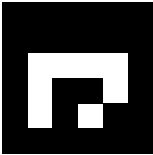


FFPE / BLOCK PHOTOGRAPHY

CAPTURE MAT v1 | A4 | Print page 1 at 100% / Actual size

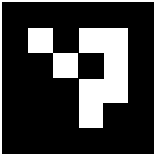

TOP / LABEL END


Keep all four markers clear



BLOCK ZONE / 65 x 55 mm



mm


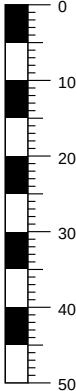


Place tissue face up; label end toward TOP.

BLOCK ID / One matching label per image




Use a barcode / QR label with readable ID text.

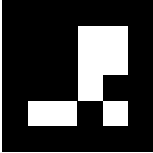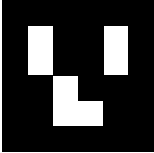

1x lens / 3:4 portrait / no digital zoom

Move phone; keep frame fully in view.

PRINT CHECK: this line must measure 50 mm




Matte paper | Flat support | Guide on page 2

Frame: 150 x 200 mm | ArUco DICT_4X4_50 | Tags: 20 mm, 1-cell border

Photograph blocks in five steps

FFPE photography mat | Quick operator guide

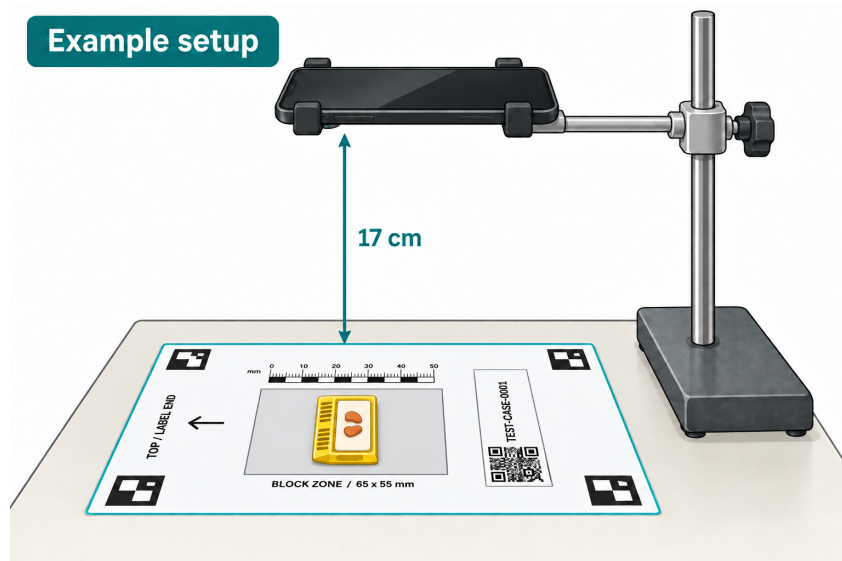

Example setup: Secure the phone screen up, with its rear lens clear and facing the mat. 17 cm is an example height; adjust to include the full mat. Illustration only.

1. Print the mat

Print page 1 on the matching paper size at **100% / Actual size**. Turn off Fit and Shrink. Check the 50 mm line with a ruler. Place the mat flat on a firm surface.

2. Match the block and label

Place one block in the gray zone, with its label end toward **TOP**. Put the matching human-readable ID and QR code in the ID box. Keep the white space around the QR code clear. Change the block and its label together.

3. Take two photos per block

Photograph the **identifier / cassette side**, then the **tissue side**. Keep the same matching ID and QR label visible in both photos. After turning the block over, keep its label end toward TOP.

4. Keep the whole mat visible

Use the rear camera at **1x**, directly above and parallel to the mat. Keep all four corner markers and the ruler in frame. Use even light, avoid glare, and tap to focus. Move the phone back if needed; do not crop, mirror or digitally zoom.

5. Transfer and ingest

Transfer the original photos to an approved folder on the laptop using your approved transfer method. Open Spatial Prep through its local server and choose **Ingest photo folder**. The app pairs by QR and aligns the photos. Review any failures before drawing scoring shapes. Use JPEG, PNG or WebP; convert HEIC locally if needed.

PHI: approved capture

Block IDs and photographs may contain **protected health information (PHI)**. Use an institution-approved camera and capture workflow, for example a managed smartphone with an approved work app that is authorized to capture and store PHI.

Having a work app installed does not make the personal Camera app, camera roll or cloud backup approved. Use only approved storage and transfer routes; avoid personal photo sync, messaging and email. Follow local retention and deletion rules for phone copies, laptop files and exported reports. If approval is unclear, use dummy material until your institution confirms the workflow.

Scale reminder: The mat measures the paper plane. The raised block face can have a different scale; use the validated capture setup and scoring tolerance for your workflow.